

Endotoxin removal by 3M™ Emphaze™ AEX Hybrid Purifier

Rituparna Sengupta, Ph.D and Rebecca Hochstein, Ph.D
3M Separation and Purification Sciences Division, 3M Center, St. Paul, MN



3M™ Emphaze™ AEX Hybrid Purifier Laboratory Capsules for Endotoxin Removal in Pre-clinical and Early Stage Recombinant Proteins

The 3M™ Emphaze™ AEX Hybrid Purifier is a synthetic, single-use anion exchange product that utilizes quaternary amine (Q) functionalized non-woven hydrogel media for clarification and purification of process fluids in biopharmaceutical purification. 3M Emphaze AEX Hybrid Purifier redefines and simplifies biologics manufacturing processes by utilizing chromatography for particle control and reduction of DNA, host cell proteins, endotoxin, and bioburden. The anion exchange chromatography of the functional nonwoven in 3M Emphaze AEX Hybrid Purifier provides substantial reduction of soluble impurities and of insoluble particles below 0.1 microns for superior turbidity reduction.

The introduction of two small laboratory capsule sizes (BV1R and BV0.3R) for the 3M Emphaze AEX Hybrid Purifier enables the use of this device for endotoxin removal in laboratory applications using peristaltic pumps, as well as using syringe filtration. The laboratory capsules allow for fast and easy processing of small volumes of early stage



drug candidates during pre-clinical and early stage clinical studies in animal-based systems, where endotoxin removal is a requirement in order to avoid eliciting an unwanted immune response while still being able to observe drug efficacy.

Application Data for Endotoxin Removal from Protein Solutions Using 3M™ Emphaze™ AEX Hybrid Purifier

Protein solutions in phosphate-buffered saline (PBS) buffer were used to represent early stage drug candidates in pre-clinical trials. Bovine serum albumin (BSA, pI = 4.5) and a biosimilar IgG1 (pI = 8.5) were chosen as model proteins to demonstrate endotoxin removal by 3M Emphaze AEX Hybrid Purifier over a practical range of isoelectric points that models the diversity of biopharmaceutical drug candidates. Proteins were spiked with 10² and 10³ EU/mL of stock endotoxin from a lysed and clarified *E. coli* culture.

3M Emphaze AEX Hybrid Purifier successfully reduced endotoxin to below WFI (water for injection) levels (<0.25 EU/mL)¹ in protein solutions at challenge loads of >1 kg/L. Flow rates of 1, 5, and 10 BV/min all reduced endotoxin levels to below WFI, demonstrating that the performance was unaffected by flow rate (Fig. 1). Syringe filtration (~ 23 BV/min) using 3M Emphaze AEX Hybrid Purifier also removed endotoxin to below WFI levels (Fig. 2).

Robust endotoxin removal to below WFI with 90 – 100% protein recovery across isoelectric points and flow rates by 3M™ Emphaze™ AEX Hybrid Purifier

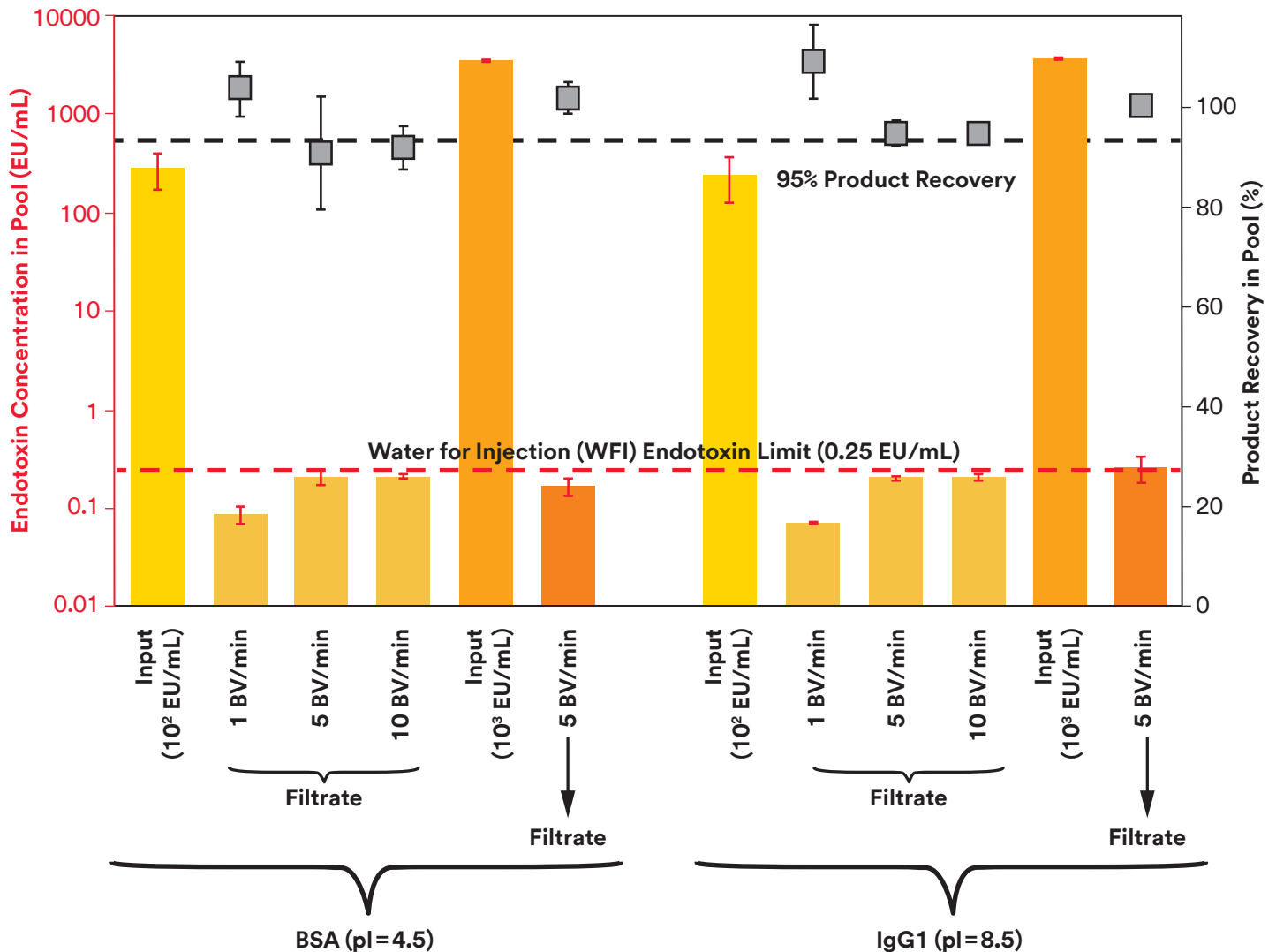


Figure 1. Endotoxin removal by 3M™ Emphaze™ AEX Hybrid Purifier for two different proteins (BSA and IgG1) at different flow rates (1 BV/min, 5 BV/min and 10 BV/min) using two different input concentrations for endotoxin (10² and 10³ EU/mL). Protein isoelectric points (pIs) were experimentally determined by isoelectric focusing. In all instances, endotoxin was reduced to below WFI (<0.25 EU/mL) levels and protein recovery was greater than 95%.

Endotoxin removal by 3M™ Emphaze™ AEX Hybrid Purifier using syringe filtration

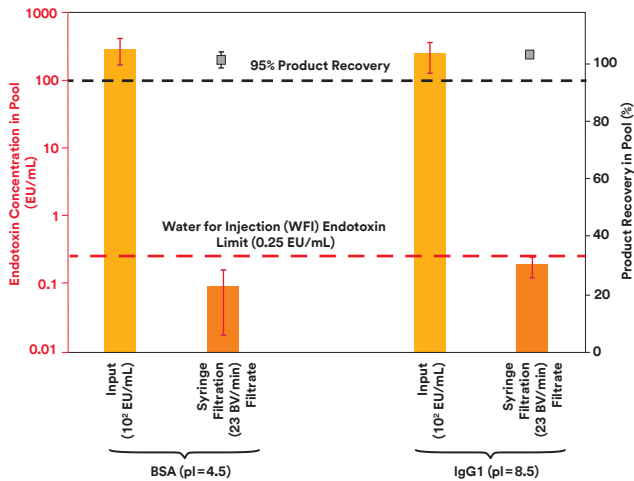


Figure 2. Syringe filtration using 3M™ Emphaze™ AEX Hybrid Purifier demonstrates endotoxin reduction to below WFI levels for an endotoxin input of 10² EU/mL at protein challenge loads of greater than 1 kg/L.

Comparison of 3M™ Emphaze™ AEX Hybrid Purifier with a Standard AEX Pre-packed Resin Column with Q-chemistry

Endotoxin removal was performed using a standard quaternary amine pre-packed AEX resin column as a comparison to 3M Emphaze AEX Hybrid Purifier. The column showed decreased endotoxin reduction at flow rates higher than 1 BV/min (Fig. 3) and levels above WFI.

3M Emphaze AEX Hybrid Purifier is ideal for endotoxin removal applications in a laboratory setting, especially for high throughput processing of protein samples for screening

during pre-clinical and early-stage clinical testing. The insensitivity to flow rate, unlike a diffusion-limited resin column, allows for the use of less complex pumping systems, such as a peristaltic pump or syringe filtration, enabling multiple samples to be processed simultaneously or in quick succession. The high capacity of 3M Emphaze AEX Hybrid Purifier allows for high protein loads, and scalable devices allow for easy transition to pilot or production scale processes. 3M Emphaze AEX Hybrid Purifier is a robust and simple solution to the challenging problem of endotoxin removal to help accelerate the preclinical development of biopharmaceutical drug candidates.

Endotoxin removal performance is reduced at high flow rates when using a standard AEX pre-packed resin column with Q-chemistry

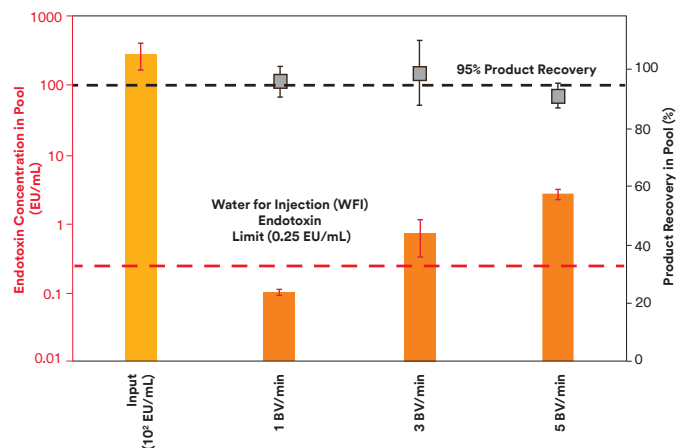


Figure 3. Endotoxin removal by a standard AEX pre-packed resin column with Q-chemistry, using BSA and an input endotoxin concentration of 10² EU/mL, shows endotoxin reduction to below WFI only when the column is operated at a flow rate of 1 BV/min. At higher flow rates, endotoxin levels were above WFI levels.

References

- 1 United States Pharmacopeia and National Formulary.
USP 29-NF 24 chapter 1231. Water for Pharmaceutical Purposes.



Intended Use: 3M™ Emphaze™ AEX Hybrid Purifier single-use filter products are intended for use in biopharmaceutical processing applications of aqueous based pharmaceuticals (drugs) and vaccines in accordance with the product instructions and specifications, and cGMP requirements, where applicable.

Restrictions on Use: 3M advises against the use of these 3M products in any application other than the stated intended use(s), since other applications have not been evaluated by 3M and may result in an unsafe or unintended condition. Do not use in any manner whereby the 3M product, or any leachable from the 3M product, may become part of or remains in a medical device that is regulated by any agency, and/or globally exemplary agencies, including but not limited to: a) FDA, b) European Medical Device Directive (MDD), c) Japan Pharmaceuticals and Medical Devices Agency (PMDA) or in applications involving permanent implantation into the body; Life-sustaining medical applications; Applications requiring food contact compliance.

Technical Information: The technical information, guidance, and other statements contained in this document or otherwise provided by 3M are based upon records, tests, or experience that 3M believes to be reliable, but the accuracy, completeness, and representative nature of such information is not guaranteed. Such information is intended for people with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. No license under any 3M or third party intellectual property rights is granted or implied with this information.

Product Selection and Use: Many factors beyond 3M's control and uniquely within user's knowledge and control can affect the use and performance of a 3M product in a particular application. As a result, end-user is solely responsible for evaluating the product and determining whether it is appropriate and suitable for end-user's application, including completing a risk assessment that considers the product leachable characteristics and its impact on drug safety conducting a workplace hazard assessment and reviewing all applicable regulations and standards (e.g., OSHA, ANSI, etc.). Failure to properly evaluate, select, and use a 3M product and appropriate safety products, or to meet all applicable safety regulations, may result in injury, sickness, death, and/or harm to property.

Warranty, Limited Remedy, and Disclaimer: Unless a different warranty is specifically stated on the applicable 3M product packaging or product literature (in which case such warranty governs), 3M warrants that each 3M product meets the applicable 3M product specification at the time 3M ships the product. 3M MAKES NO OTHER WARRANTIES OR CONDITIONS, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OR CONDITION OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR ARISING OUT OF A COURSE OF DEALING, CUSTOM, OR USAGE OF TRADE. If a 3M product does not conform to this warranty, then the sole and exclusive remedy is, at 3M's option, replacement of the 3M product or refund of the purchase price.

Limitation of Liability: Except for the limited remedy stated above, and except to the extent prohibited by law, 3M will not be liable for any loss or damage arising from or related to the 3M product, whether direct, indirect, special, incidental, or consequential (including, but not limited to, lost profits or business opportunity), regardless of the legal or equitable theory asserted, including, but not limited to, warranty, contract, negligence, or strict liability.



3M Purification Inc.
3M Separation and Purification Sciences Division
400 Research Parkway, Meriden, CT 06450 USA

Phone 1-800-243-6894 1-203-237-5541

Web 3M.com/bioprocessing

3M and Emphaze are trademarks of 3M Company. All other trademarks are property of their respective owners.

Please recycle. Printed in USA © 3M 2020.
All rights reserved.