



**3M™ Attest™ Vaporized
Hydrogen Peroxide Tri-Metric
Chemical Indicator 1348/1348E**

Monitor three key variables with one simple innovation.

The first VH2O2 CI to independently monitor
all three variables – time, temperature,
and VH2O2 concentration.

For the first time, you can get a higher level of assurance with VH2O2 monitoring as you currently expect with steam monitoring. The innovative 3M™ Attest™ Tri-Metric Chemical Indicator, Type 4, streamlines your workflow by displaying clear sterilization results for heat and moisture sensitive instruments.



**Detects failed VH2O2
sterilization cycles** more accurately than any other CI tested.*1



Easy-to-read 'accept' and 'reject' zones. Clearly changes from blue to pink.²



The most colour stable of all VH2O2 CIs tested* – won't change colour before sterilization under ambient light.³



Proprietary protective film prevents ink smear on surgical instruments.⁴



Leading the way for better patient care


The first FDA cleared CI for VH2O2 sterilization categorised by the FDA as a multivariable CI for chemical vapour sterilization that independently monitors three critical variables.⁵

- time**
- temperature**
- VH2O2 concentration**


Expect a higher level of performance.

VH2O2 sterilization is designed to perform at lower temperatures to process heat and moisture sensitive instruments and devices. With the 3M™ Attest™ VH2O2 Tri-Metric Chemical Indicator, Type 4, you get sterilization monitoring results you can count on.

- ▶ Easy-to-find among instruments in a pack
- ▶ Clear results at a glance. Moving front technology eliminates time spent interpreting colour changes.
- ▶ 2D barcode contains lot numbers and expiration dates to scan to your hospital’s current instrument tracking system
- ▶ Blue to pink colour change is consistent on all 3M chemical indicators for vaporized hydrogen peroxide

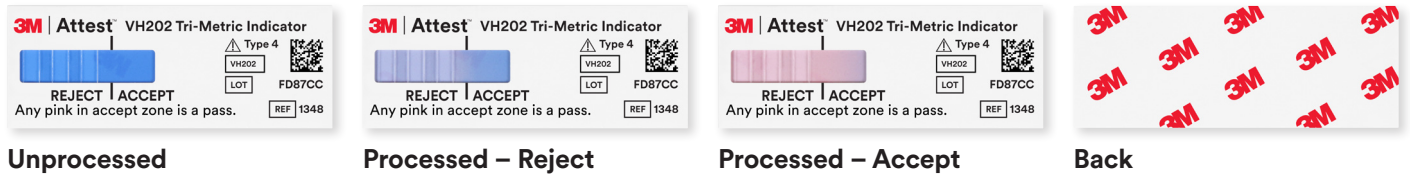

Product	Product name	Size	Items/bag	Per case
1348	3M™ Attest™ Vaporized Hydrogen Peroxide Tri-Metric Chemical Indicator	5.08 cm	250	4 bags
1348E	3M™ Attest™ Vaporized Hydrogen Peroxide Tri-Metric Chemical Indicator	10.16 cm	250	4 bags
1295	3M™ Attest™ Rapid Readout Biological Indicator (pink cap)	–	30	4 bags
1228	3M™ Comply™ Hydrogen Peroxide Indicator Tape	1.90 cm	–	24 rolls

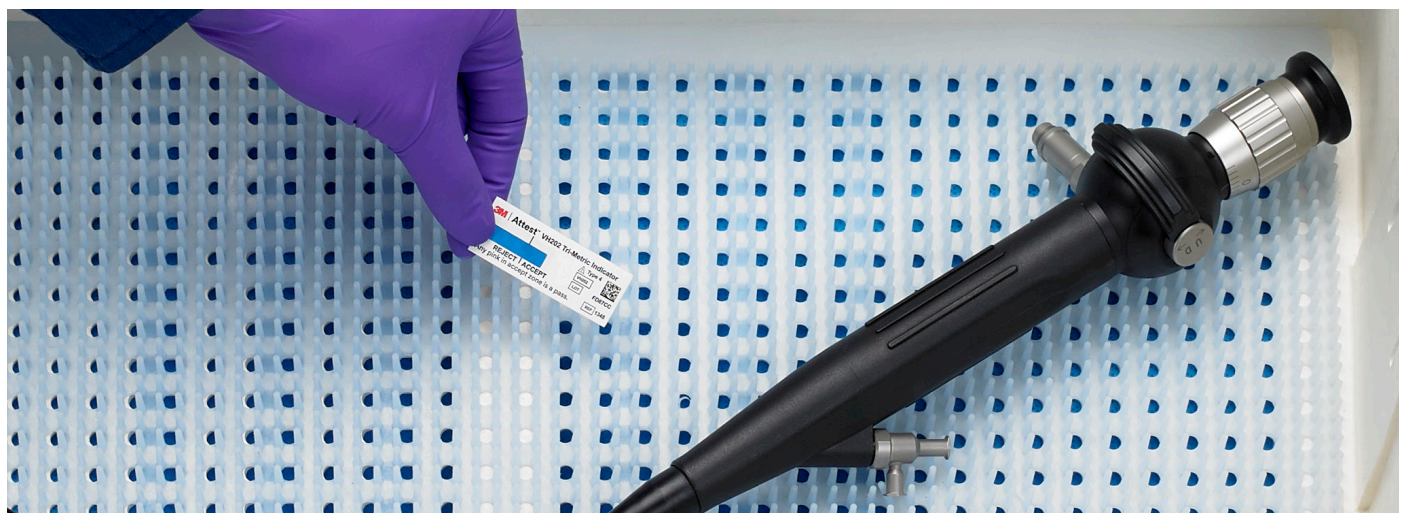

For more information, or to request a sample, contact your sales representative or visit go.3M.com/TriMetric

3M United Kingdom PLC
 Charnwood Campus
 10 Bakewell Road
 Loughborough
 LE11 5RB
 01509 611611

3M Ireland
 The Iveagh Building
 Carrickmines Park
 Carrickmines
 Dublin 18
 00 353 (01) 280 3555

1 EM-05-692856, EM-05-692858
 2 EM-05-424823
 3 EM-05-692861
 4 EM-05-638321
 5 EM-05-673143

3M, Attest and Comply are trademarks of 3M Company. © 3M 2020. OMG57577.

